
SonicBrush Clean30

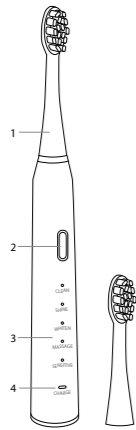EN - INSTRUCTIONS


PACKAGE CONTENTS

- Sonic brush TrueLife SonicBrush Clean30
- 2x Standard heads
- USB-C cable
- Instructions
- Safety instructions

PRODUCT DESCRIPTION

1. Interchangeable head
2. ON / OFF button and mode change
3. Modes indicator
4. Charging indicator

TOOTHBRUSHING

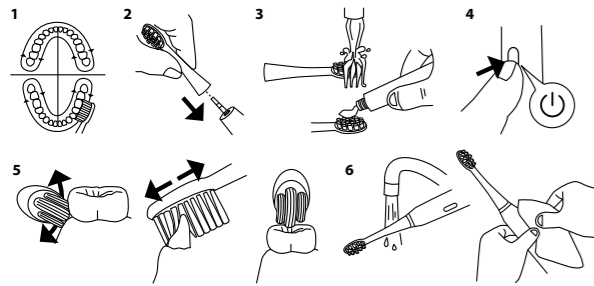1. For maximum efficiency, divide the cleaning into 4 sections and devote 30 seconds to each of them. After 30 seconds, the brush will take a short pause and alert you to move to the next section. After 2 minutes, the brush turns itself off. Cleanse at least twice a day, in the morning and in the evening.
2. Fit the brush head firmly onto the metal shaft with the bristles facing the front of the handle.
3. Moisten the head, apply toothpaste on it and rub it on the teeth.
4. Switch on the brush with the ON / OFF button. To select another mode, press the button again and select one of the five modes. If you would like to end cleaning earlier, hold the ON/OFF button.
5. In each section, start by cleaning the outside of the teeth, then the inside of the teeth, and finally clean the surfaces. Hold the brush at an angle so that the bristles gently touch the teeth and the edges of the gums at the same time. When cleaning, do not press the brush and do not make circular movements like with a manual brush, but move slowly from one tooth to another.
6. Rinse the brush head and body under running water and then wipe them with a cloth or towel.

CLEANING MODES

- CLEAN – For routine cleaning and removal of dental plaque
- SHINE – To remove coffee and tea stains and achieve a brilliant shine on teeth
- WHITEN – To remove stubborn dental plaque and whiten teeth
- MASSAGE – To stimulate the gums and improve their health
- SENSITIVE – For cleaning sensitive teeth and gums

REPLACEMENT OF HEADS

We recommend changing the heads every 3 months. After this period of use, wear occurs when the frayed and bent fibers have a reduced ability to remove dental plaque and bacteria accumulate on them.

CHARGING

1. When the brush is almost empty, the charging indicator will flash red.
2. Connect the connector of the USB cable to the socket of the brush, the other end to the adapter and the adapter to the socket. You can also charge the brush from a computer or power bank. Charge the sonic brush using the included USB-C cable. Do not use fast chargers for charging.
3. During charging, the indicator flashes green. Once the brush is charged, the charging indicator lights up green. Charging time is 3 hours.

SPECIFICATIONS

Battery life	30 days (2x a day 2 minutes)
--------------	------------------------------

Charging	USB-C kabel
Number of movements	70 000 movements/min
Modes	Clean, Shine, Whiten, Massage, Sensitive
Timer	2 min (30 s interval)
Memory	Yes
Additional heads	Clean-series
Brush dimensions	26 x 26 x 23,5 cm
Brush weight	108 g
Waterproof	IPX7
Input voltage	5 V – 1 A
Battery	600 mAh, Li-ion
Charging time	3 hours

MANUFACTURER'S RESPONSIBILITY

TrueLife products are warranted against manufacturing defects. In the event of such defects under warranty, elem6 s.r.o. shall repair or replace the defective part or product with a comparable part or product at its option. Rebuilding, altering or otherwise tampering with the product shall be deemed to exceed the scope of this warranty. This is the case even if the product is lost or damaged as a result of a defect, negligence or other error in the product. This warranty does not cover damage resulting from use, accident or normal wear and tear. elem6 s.r.o. assumes no liability for any accident, injury, death, loss or other claim related to or arising out of the use of this product. In no event shall elem6 s.r.o. be liable for incidental or consequential damages related to or arising out of the use of this product or any of its components. The return or replacement of parts and/or products may be subject to shipping, handling, replacement and/or restocking charges.

SAFETY WHEN WORKING WITH THE DEVICE

1. Read the safety information leaflet carefully.
2. Read the current version of the manual, which is available for download at: www.truelife.eu.
3. Printing errors and changes to information are subject to change.

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