LANAFORM

New Baby Scale

HOLD

TARE

PRECISION

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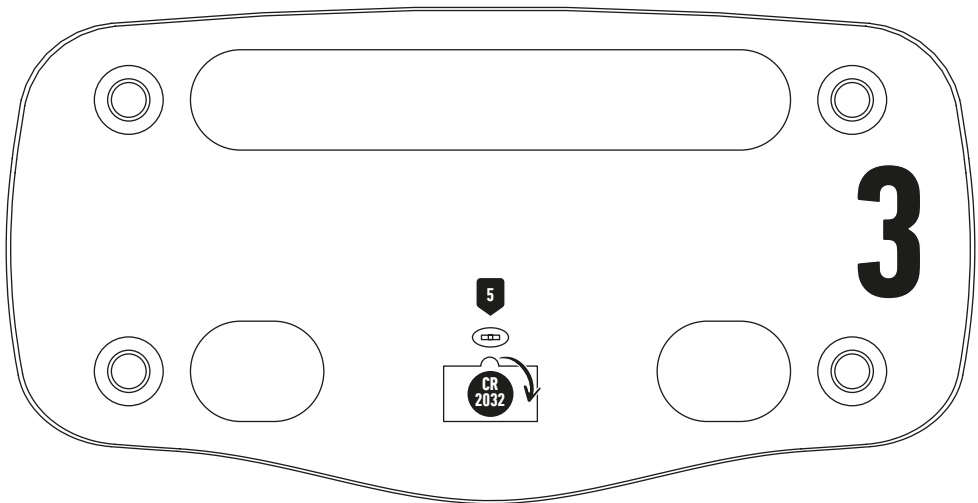
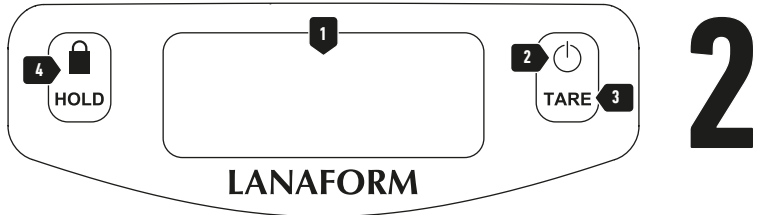
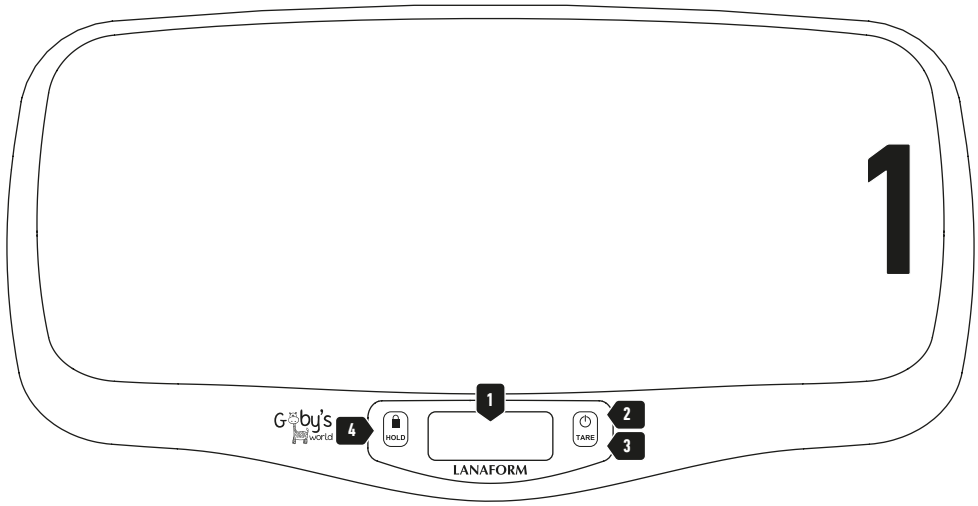
BABYWAAGE

BILANCIA PER BÉBÉS

ENGLISH

Thank you for choosing the LANAFORM New Baby Scale.

The accurate and easy to use digital New Baby Scale allows you to weigh your baby to a precision of 5 grams.

The New Baby Scale is a useful instrument for you and your baby during the first few months of their life.

IMPORTANT SAFETY MEASURES

▲ Please read all the instructions before using your product, particularly these basic safety instructions.

📷 The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

Never leave the baby unsupervised when using the New Baby Scale.

The scale is not suitable for babies weighing over 20kg. Overloading the scale may cause damage and injury.

Always place the New Baby Scale on a hard, flat, horizontal surface, away from sources of heat such as stoves, radiators etc.

Do not use this device in a room in which aerosols (sprays) are used or in a room in which oxygen is administered.

Do not use the device near wet areas, such as the bath, sink etc.

Do not immerse the device in water or any other liquid.

When using the scale, keep them away from sources of electromagnetic waves such as microwave ovens, mobile telephones and wireless networks, as this may affect the measurement.

To avoid electrocution, do not try to repair the device yourself.

In case of electrostatic failure, install the batteries again to resume normal performance.

If the screen is showing an "Err" error message, switch off the scale. Then power up again and follow the instructions.

Never leave anything on the scale platform when not in use.

When not in use, the scale switch off automatically after a few minutes. To switch off the scale more quickly, press the **ON/OFF** button for half a second.

This equipment is not intended for use by people, including children, with limited physical, sensory or mental abilities, or by people lacking experience or knowledge, unless

they are supervised or have been instructed beforehand on how to use the equipment by someone responsible for their safety. Children should be supervised to make sure they do not play with the device.

Only use the New Baby Scale for its intended purpose, as described in the manual.

COMPOSITION OF THE PRODUCT [📄 1, 2 & 3]

- 1 Screen
- 2 **ON/OFF** button
- 3 **TARE** (Deduction of initial weight)
- 4 **HOLD** (Determination of weight at a specific moment)
- 5 **KG/LB/OZ** (Measurement unit)

TECHNICAL FEATURES

- Capacity: 0.05 to 18 kg.
- Scale markings: 5 grams.
- Tare weight function.
- HOLD function.
- LCD screen
- Change measuring unit function (kg/lb/oz)
- Low battery alert: the "Lo" icon appears when the battery is low.
- Automatic switch-off.
- Batteries: 1 x CR2032 (supplied).
- Screen size: 72 (L) x 28 (W) mm.
- Device size: 560 (L) x 330 (W) x 53 (H) mm.
- Net weight of device: 1.56 kg.
- Storage environment: -20°C to 60°C (humidity: 25%-85%).
- Operating environment: 0°C to 40°C (humidity: 20%-95%).
- Measurement precision:
 - 0.05 kg – 10 kg (±0,005 kg).
 - 10 kg – 15 kg (±0,010 kg).
 - 15 kg – 20 kg (±0,020 kg).

PREPARATION AND BATTERY REPLACEMENT [📄 3]

Installing the battery

- 1 Remove the battery compartment cover underneath the scale.
- 2 Remove the plastic protection from the CR2032 battery or insert a new CR2032 battery if necessary.
- 3 Replace the cover.
- 4 Place the scale on a hard, flat surface before weighing your baby.
- 5 Before placing the baby on the scale, make sure that the surface of the scale is completely dry and not slippery. You can also place a towel on the scale.

"KG/LB/OZ" UNIT SELECTOR FUNCTION [📄 3]

- 1 Before placing the baby on the scale, press this button to select the measuring unit: **KG/LB/OZ**.

- 2 The measuring unit will be displayed the next time the scale is switched on.

NORMAL WEIGHING FUNCTION

- 1 Press the **ON/OFF** button to switch on the scale.
- 2 **0.000 kg (lb or oz (pounds: ounces))** depending on the measuring unit selected) will be displayed on the screen.
- 3 Place the baby on the scale. If the baby moves too much, select the **HOLD** function which will set the baby's weight when they are more still. When the "Hold" icon stops flashing, the weight will remain the same for a further 20 seconds.
- 4 "Err" will be displayed if the scale is overloaded.

WEIGHING FUNCTION WITH A TOWEL

- 1 Press the **ON/OFF** button to switch on the scale.
- 2 **0.000 kg (lb or oz (pounds: ounces))** depending on the measuring unit selected) will be displayed on the screen.
- 3 Place a towel on the scale. The weight of the towel will be displayed.
- 4 When the "Hold" icon appears and stops flashing, press the **TARE** button.
- 5 The "Tare" icon and **0.000 kg (lb or oz (pounds: ounces))** depending on the measuring unit selected) will be displayed on the screen again.

🕒 Leave the towel on the scale. If you take it off the scale, the screen will show a negative weight.

- 6 Place the naked baby on the scale. Repeat the operation with the **HOLD** button as explained above in the "normal weighing" section, if necessary.

MAINTENANCE AND STORAGE

Use a soft, damp cloth to clean the New Baby Scale.

Do not use abrasive detergents.

In case of prolonged lack of use of the device, remove the batteries.

Store the device in a dry, safe place.

🗑️ ADVICE ON WASTE DISPOSAL

The packaging consists entirely of environmentally-friendly materials that can be taken to your local waste sorting centre to be used as secondary materials. Cardboard can be placed in a paper collection container. Packing film must be taken to your local sorting and recycling centre.

When you are no longer using the device, dispose of it in an environmentally responsible manner and in compliance with the statutory requirements.

LIMITED WARRANTY

LANAFORM guarantees that this product will be free from any material and workmanship defects for a period of two years from the date it was purchased, other than as specified below.

The LANAFORM warranty does not cover damage caused by normal wear and tear on this product. In addition, the warranty on this LANAFORM product does not cover damage caused as a result of any abusive, inappropriate or incorrect use, accident, attachment of any unauthorised accessory, modification of the product or any other situation of whatever kind, beyond the control of LANAFORM.

LANAFORM shall not be held liable for any type of incidental, consequential or special damage.

Any implied warranties of the product's fitness for purpose are limited to a period of two years from the original purchase date on presentation of a copy of the proof of purchase.

Upon receipt, LANAFORM will repair or replace your device, as appropriate, and return it to you. The warranty may only be used at the LANAFORM Service Centre. Any maintenance and repair work carried out on this product by any person other than the LANAFORM Service Centre will invalidate this warranty.

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